



Belegung

FS-Innopeg-1000

ID 0,25mm, df 0,25µm, Länge 30m

Max. Temperatur: 240/260°C

Dimensionen

30 m x 0,25 mm ID, df 0,15 µm

Artikel Nr.

[22516730](#)

Probe

Hexanderivate

Substanzen: Hexanal / Cyclohexylamin / Cyclohexanon / 1-Hexanol / Cyclohexanol

Ofen


80°C isotherm

Trägergas

0,6 bar N₂

Detektor

FID 340°C

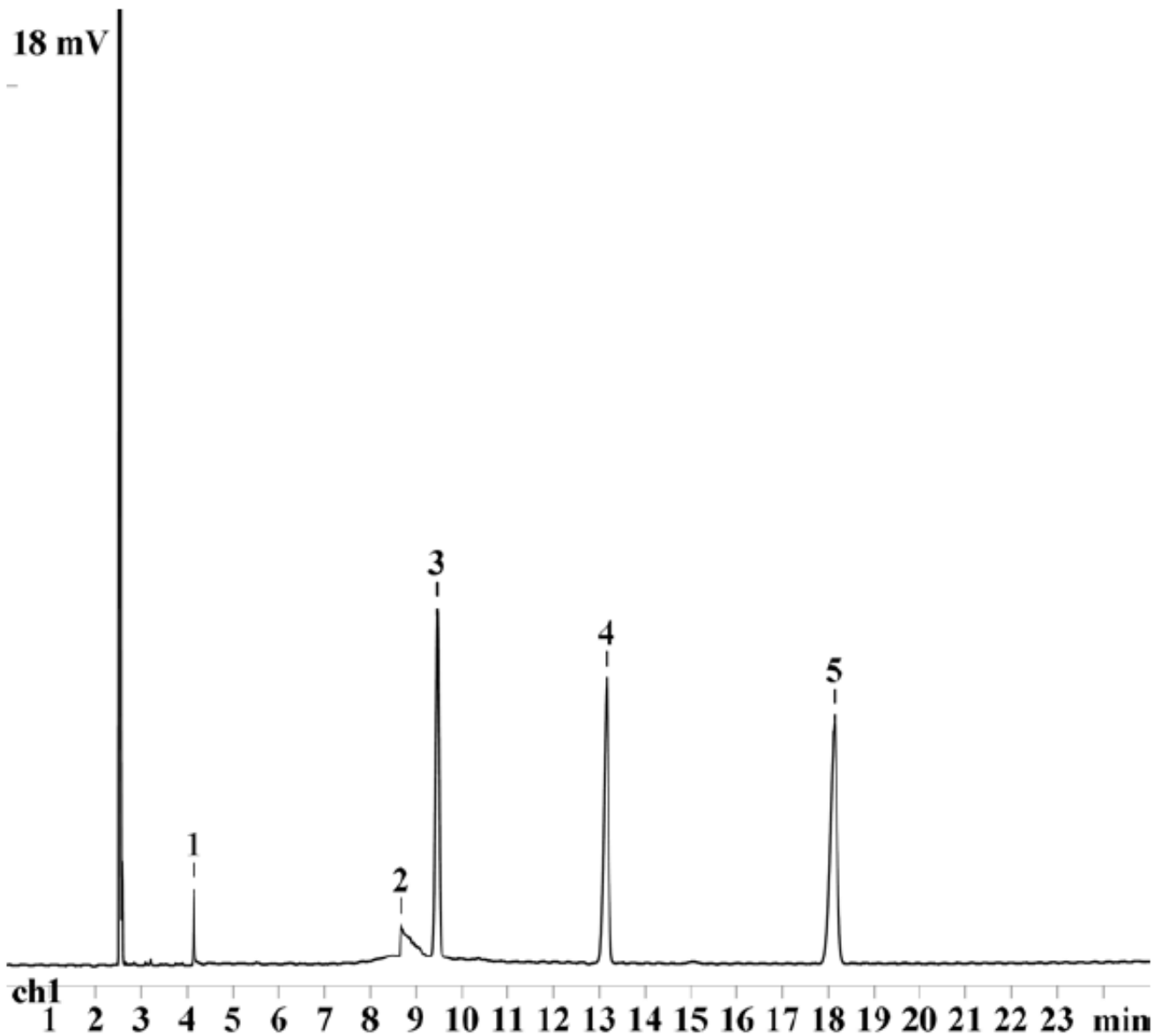

1.) Hexanal

4.) 1-Hexanol



2.) Cyclohexylamin

5.) Cyclohexanol

3.) Cyclohexanon