



Belegung

FS-Innopeg-1000

ID 0,25mm, df 0,25µm, Länge 50m

Max. Temperatur: 240/260°C

Dimensionen

50 m x 0,25 mm ID, df 0,25 µm

Artikel Nr.

[22516750](#)

Probe

2- und 3-Methylbutanol in Pentan

Ofen

60°C isotherm


Trägergas

1 bar N₂

Detektor

FID 340°C

8 mV

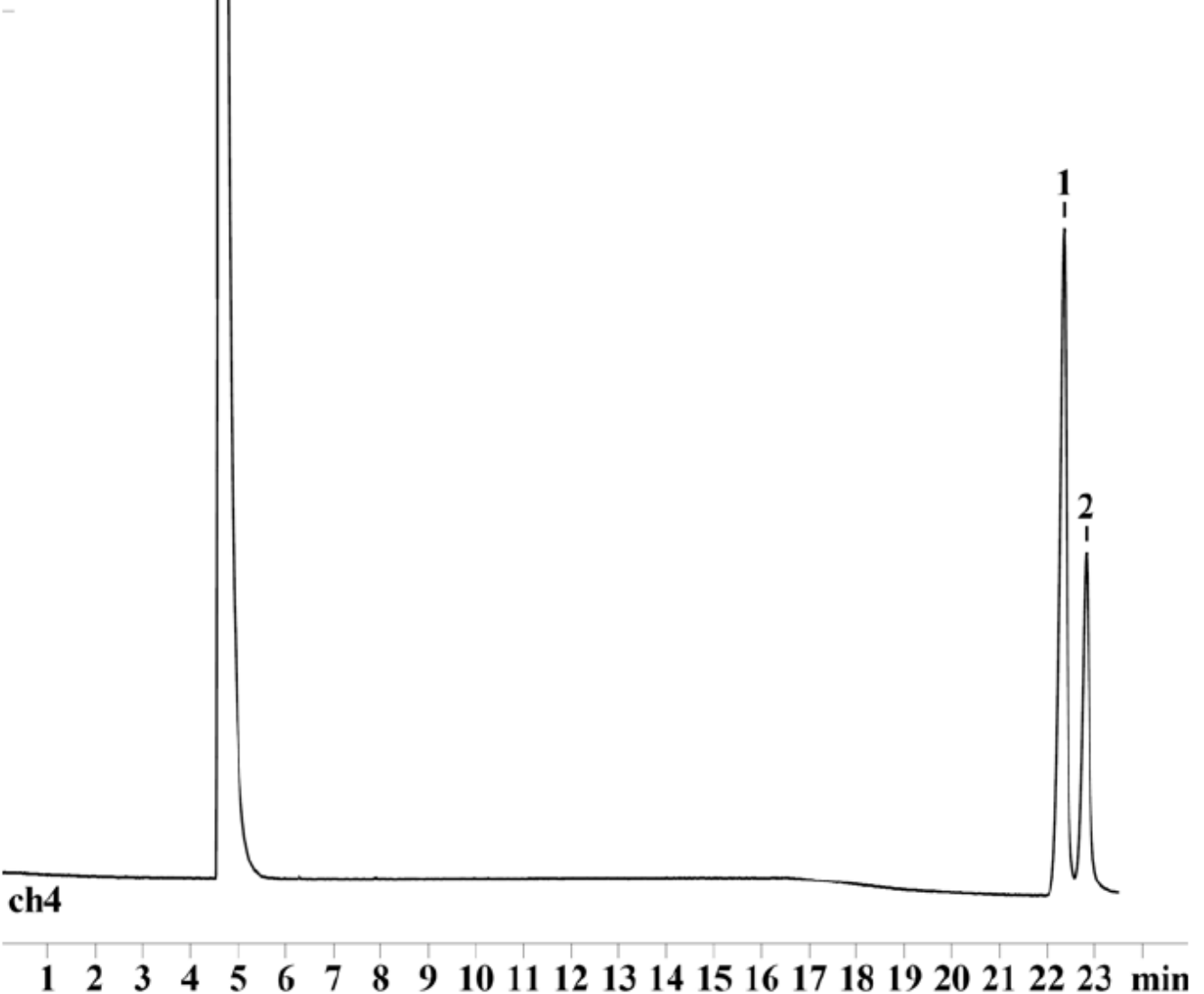

ch4

1.) 2-Methyl-1-butanol

2.) 3-Methyl-1-butanol