



Belegung

FS-Supreme-5ms

Dimensionen

60 m x 0,25 mm ID, df 1,2 µm

Artikel Nr.

[22545760](#)

Probe

Ketone in Pentan

Substanzen: Aceton / 2-Butanon / Pentanon / Isobutylmethylketon /  
Acetylaceton / Cyclohexanon

Ofen


60°C, 8°C/min, 320°C, 5min

Trägergas

1,5 bar N<sub>2</sub>

Detektor

FID 340°C

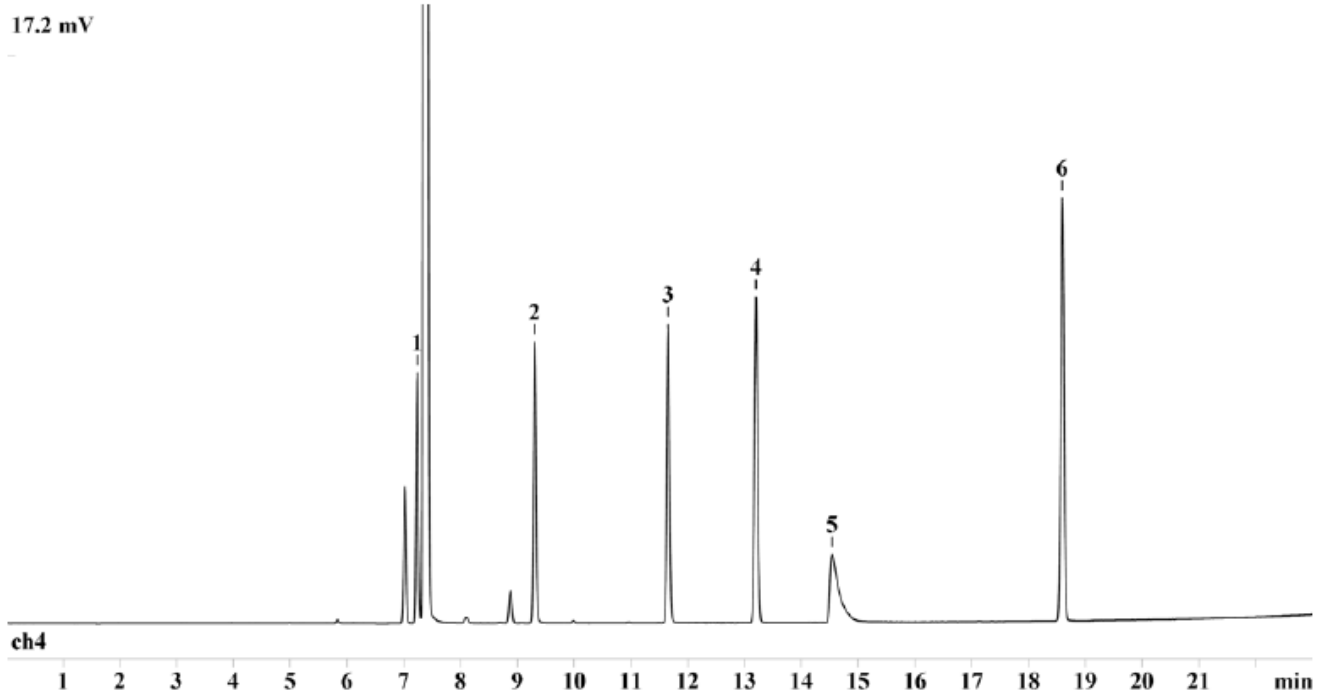

1.) Aceton

4.) Isobutylmethylketon

2.) 2-Butanon

5.) Acetylaceton

3.) Pentanon

6.) Cyclohexanon