



Belegung

FS-FS-OV-1-CB-1

Dimensionen

30 m x 0,25 mm ID, df 1 µm

Artikel Nr.

[22502330](#)

Probe

Ketone in Pentan

Substanzen: Aceton / 2-Butanon / Pentanon / Isobutylmethylketon /
Acetylaceton / Cyclohexanon

Ofen


80°C

Trägergas

0,8 bar N₂

Detektor

FID 340°C

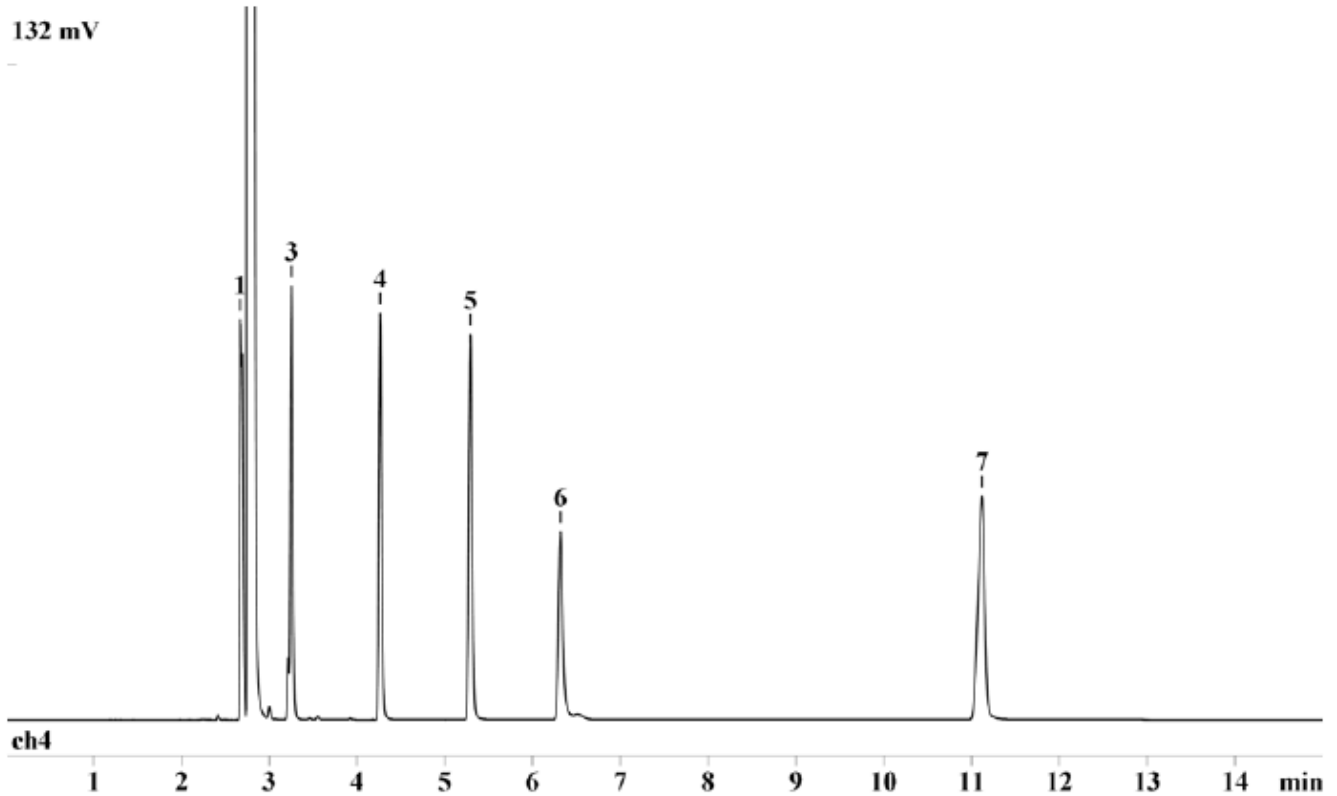

1.) Aceton

5.) Isobutylmethylketon

2.) 2-Butanon

6.) Acetylaceton

3.) Pentan

7.) Cyclohexanon

4.) Pentanon