



Belegung

FS-Innopeg-1000

ID 0,25mm, df 0,25µm, Länge 30m

Max. Temperatur: 240/260°C

Dimensionen

30 m x 0,25 mm ID, df 0,15 µm

Artikel Nr.

[22516730](#)

Probe

Aldehyde

Substanzen: Hexanal / Nonanal

Ofen


60°C, 5°C/min, 180°C, 5 min

Trärgas

0,7 bar N₂

Detektor

FID 340°C

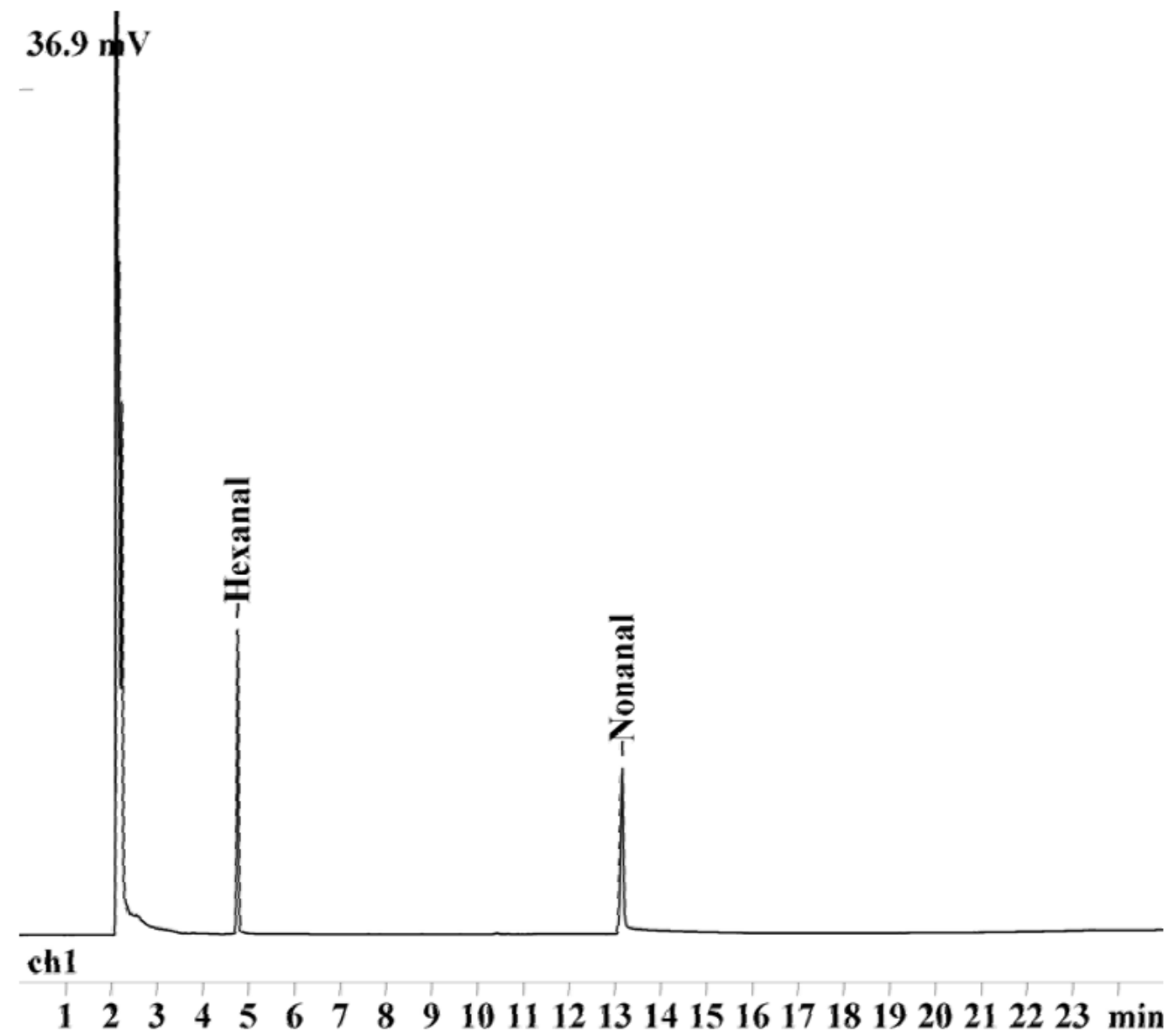

1.) Hexanal

2.) Nonanal