



Belegung

FS-Innopeg-1000

ID 0,25mm, df 0,25µm, length 50m

Max. Temperature: 240/260°C

Dimensionen

50 m x 0,25 mm ID, df 0,1 µm

Artikel Nr.

[22516750](#)

Probe

Isobutylmethylketon +

Essigsäureisobutylester in Pentan

Ofen

80°C isotherm


Trärgas

1 bar N2

Detektor

FID 340°C

336 mV

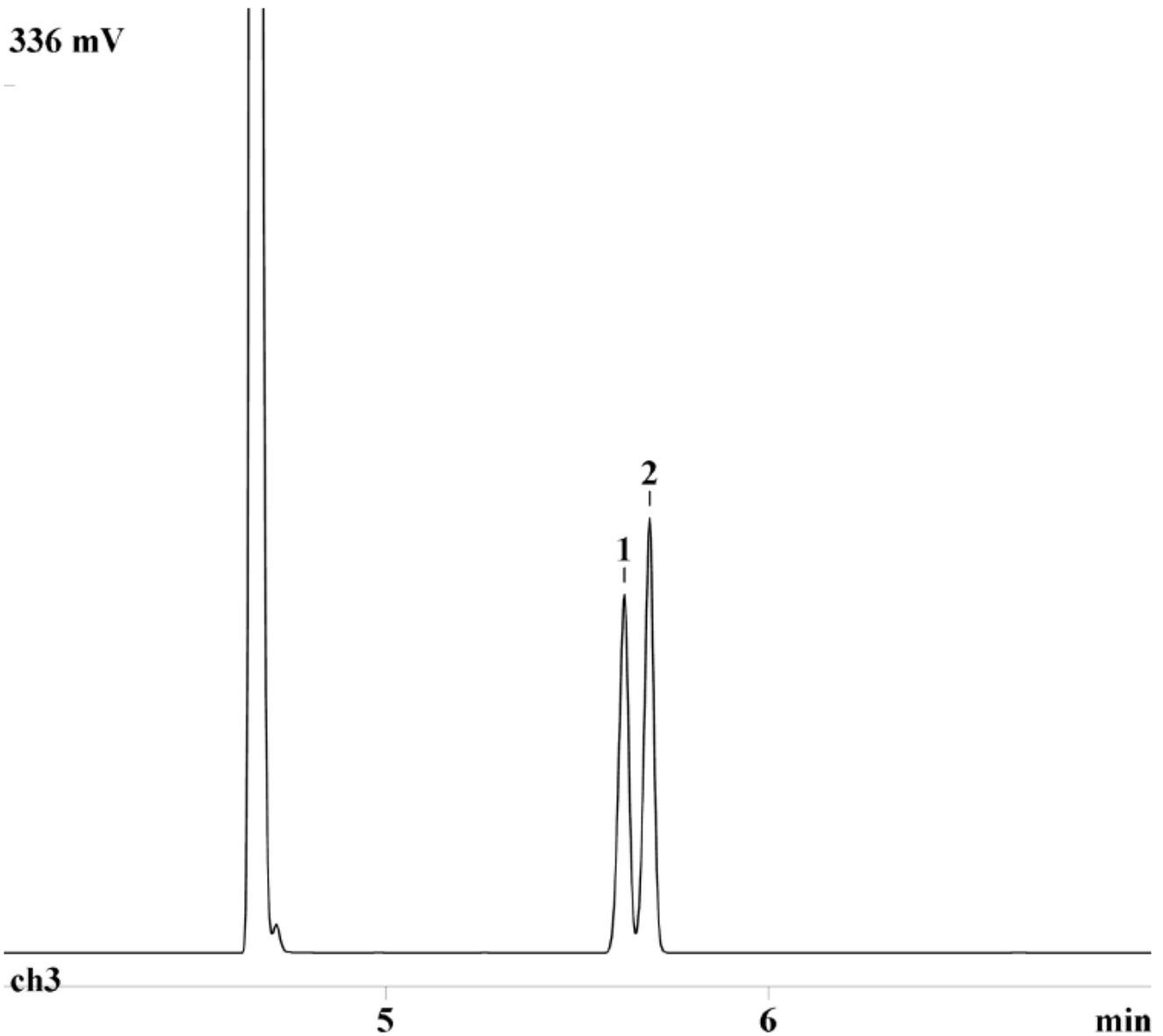

1.) Isobutylmethylketon

2.) Essigsäureisopropylester

