



Belegung

FS-Innopeg-FFAP

ID 0,25mm, df 0,25µm, Länge 30m

Max. Temperatur: 240/260°C

Dimensionen

30 m x 0,25 mm ID, df 0,25 µm

Artikel Nr.

[22517430](#)

Probe

Alkohol-Testmischung

Substanzen: n-Hexanol / n-Heptanol / n-Oktanol / n-Nonanol /
n-Decanol

Ofen


80°C; 10°C/min, 160°C, 10 min

Trägergas

1,2 bar N2

Detektor

FID 320°C

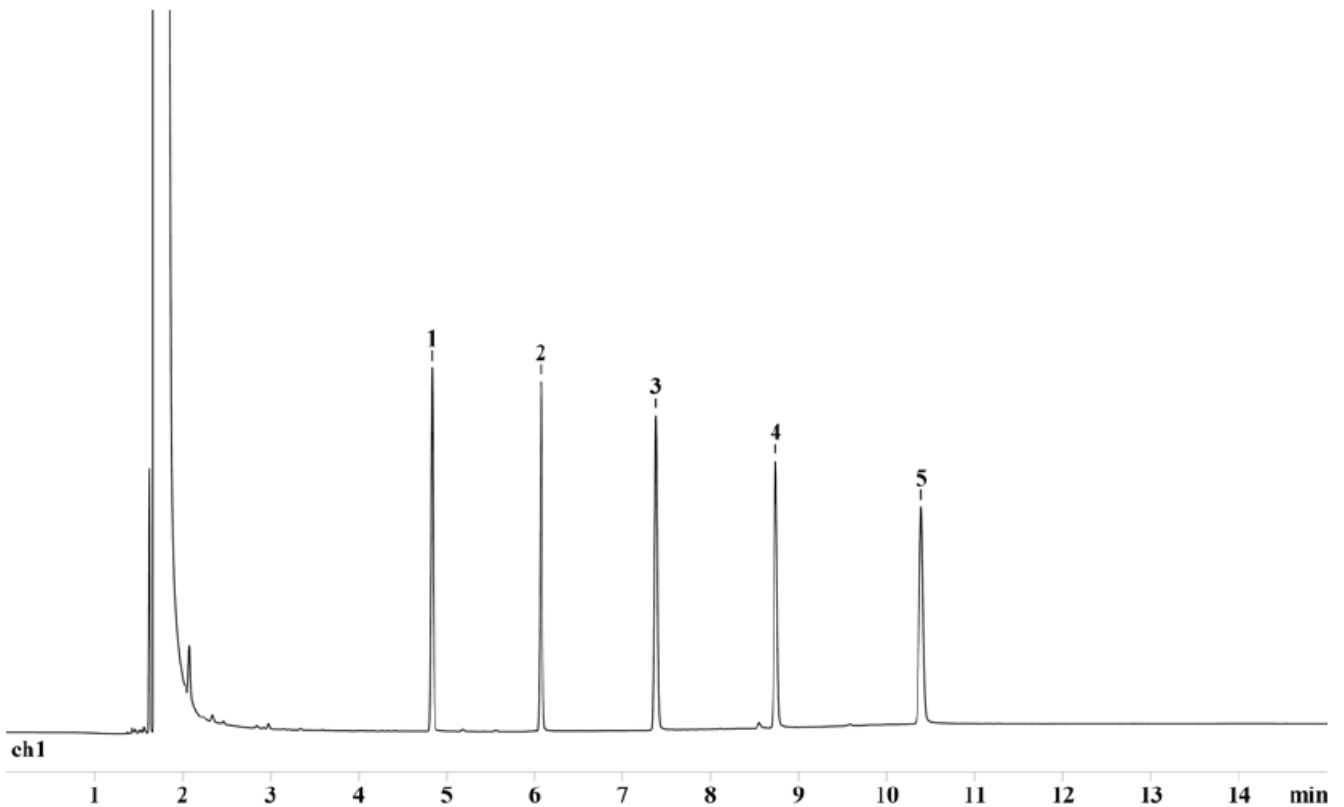

1.) n-Hexanol

4.) n-Nonanol

2.) n-Heptanol

5.) n-Decanol

3.) n-Oktanol