



Belegung

FS-OV-1 CB

Dimensionen

60 m x 0,25 mm ID, df 1  $\mu$ m

Artikel Nr.

[22502460](#)

Probe

Alkohole

Substanzen: Methanol / Ethanol / Aceton / Pentan / 1-Propanol

Ofen

45°C, 10 min, 2°C/min, 100°C


Trärgas

1 bar N<sub>2</sub>

Detektor

FID 320°C

14.5 mV

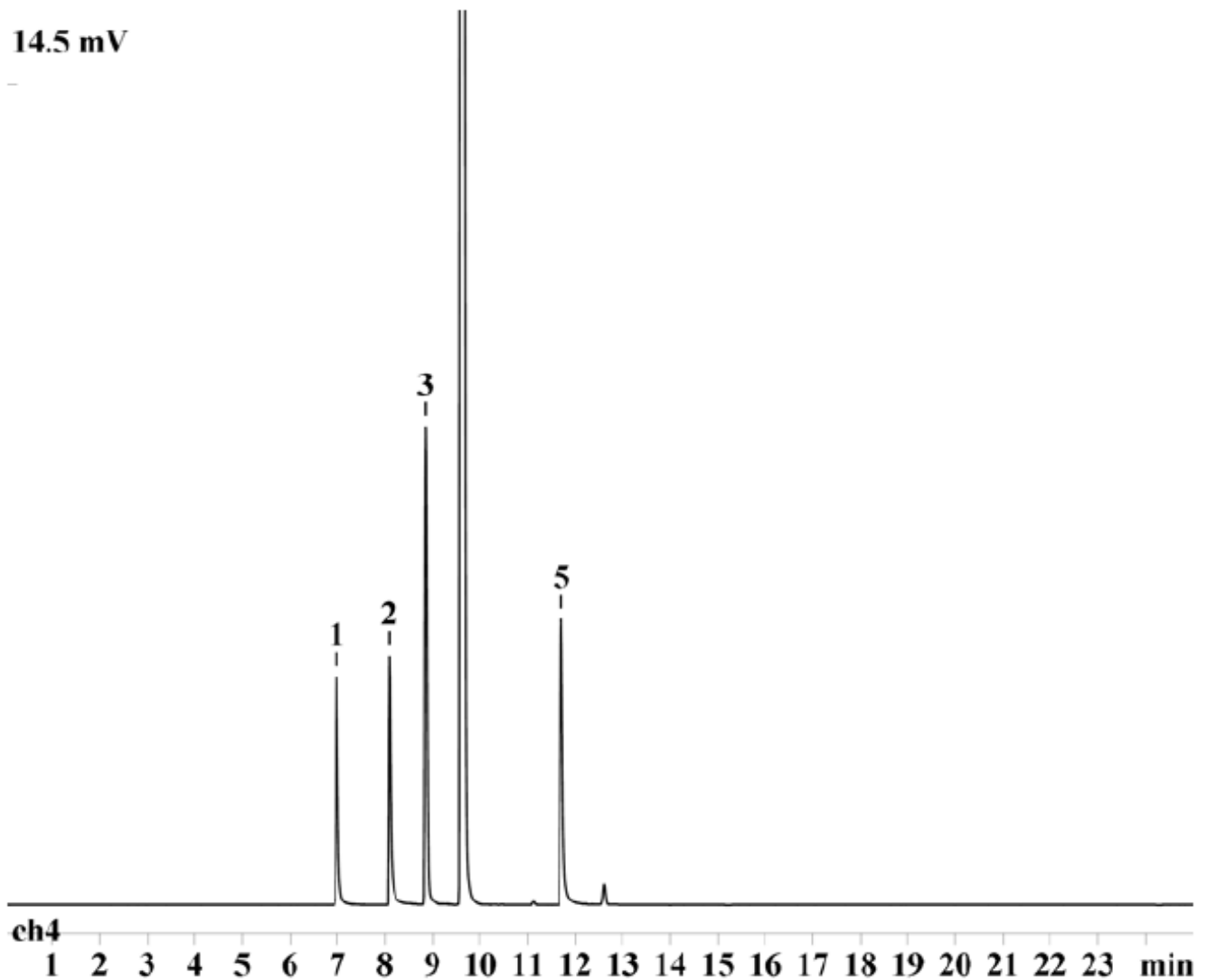

1.) Methanol

4.) Pentan

2.) Ethanol

5.) 1-Propanol

3.) Aceton