



Belegung

FS-FFAP

ID 0,25mm, df 0,25µm, length 25m

PEG 20000 modified with 2-Nitroterephthalic acid

Max. Temperature: 220°C

Dimensionen

25 m x 0,25 mm ID

Artikel Nr.

[225083](#)

Probe

Fettsäuremethylestermischung

Ofen


90-240°C mit 6°C/min

Trärgas

0,85 bar H<sub>2</sub> x Split 1:20

Detektor

FID 250°C 10 x 16

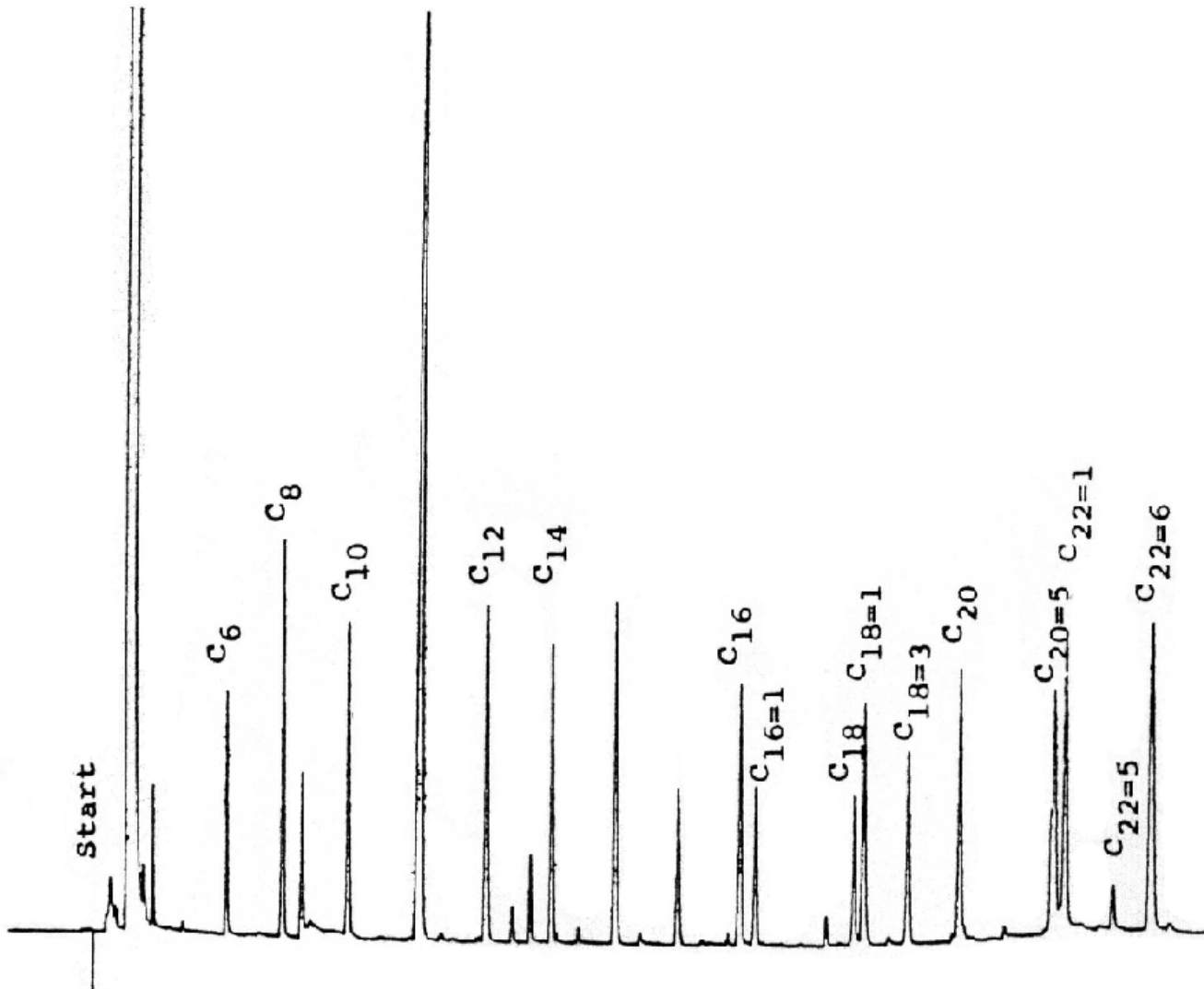

1.) FAME