

Belegung

FS-SE-54

Dimensionen

10 m x 0,32 mm ID

Artikel Nr.

[13299910](#)

Probe

Trennung von Glyceriden

Substanzen: Monoglyceride / Diglyceride / Triglyceride /
Epoxyglyceride / Dihydroxyglyceride


Ofen

160-350°C mit 7°C/min

Trägergas

Detektor

FID 380°C

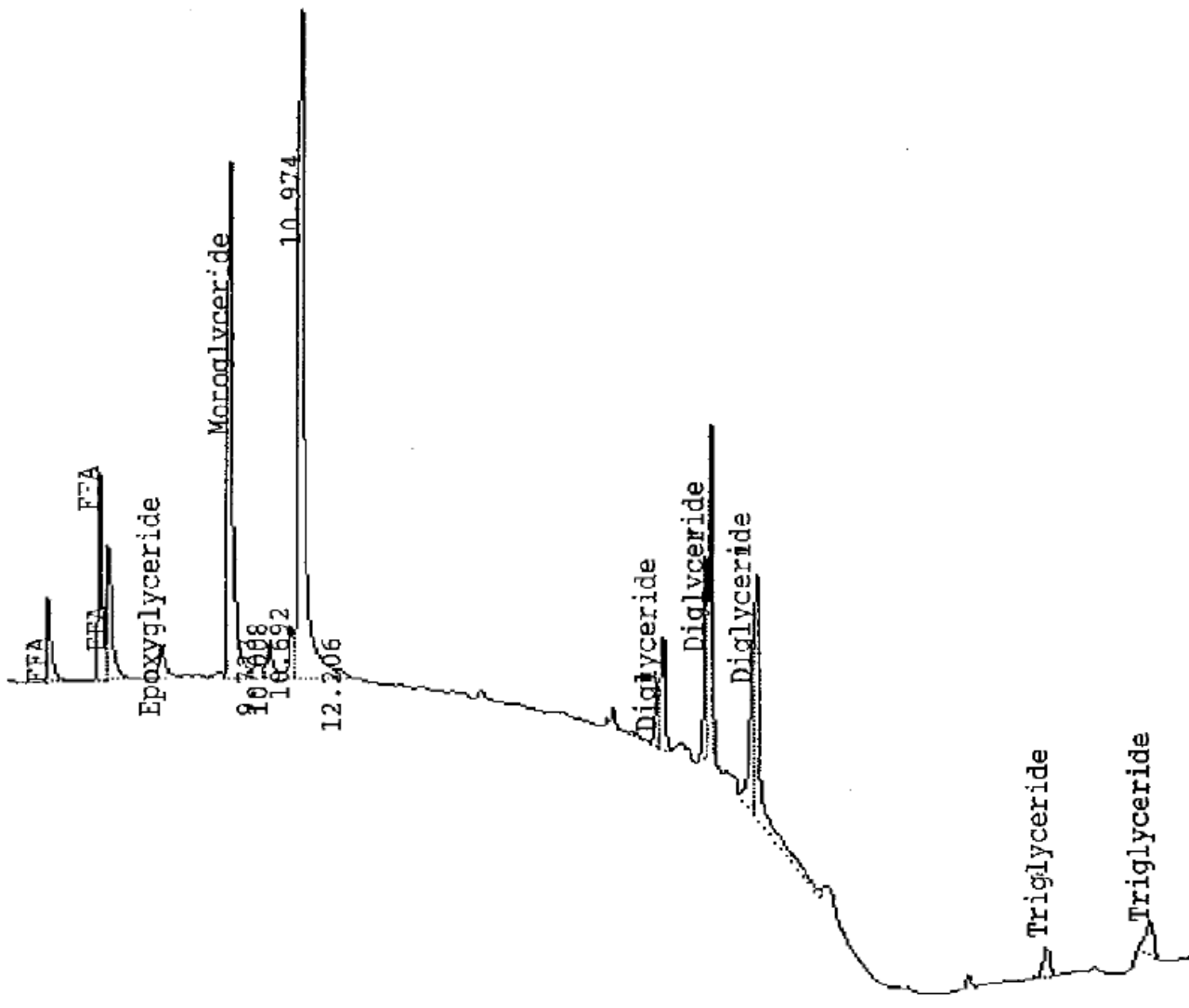

- 1.) Epoxyglyceride
- 2.) Monoglyceride
- 3.) Dihydroxyglyceride

- 4.) Diglyceride
- 5.) Triglyceride

