



Belegung

FS-SE-54

Dimensionen

2 x 50 m x 0,32 mm ID (gekoppelt)

Artikel Nr.

[232057100](#)

Probe

Trennung Alkane (C1-C6)

Substanzen: Methan / Ethan / Propan / iso-Butan / n-Butan / n-Pentan / n-Hexan

Ofen


80°C isotherm

Trärgas

1 bar N2, Split geschlossen

Detektor

FID 200°C, 1 x 16

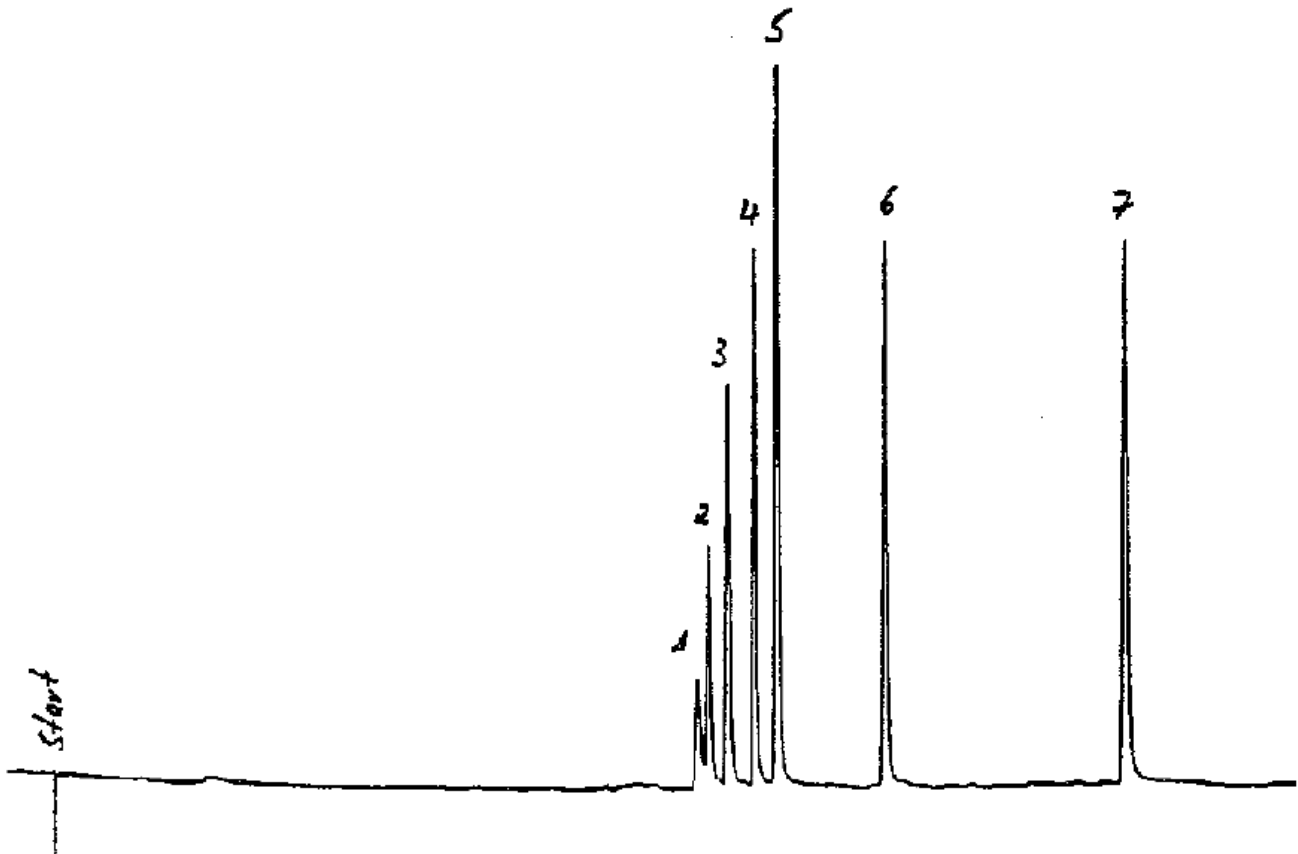

- 1.) Methan
- 2.) Ethan
- 3.) Propan
- 4.) iso-Butan

- 5.) n-Butan
- 6.) n-Pentan
- 7.) n-Hexan