



Belegung

FS-FFAP

ID 0,32mm, df 0,25µm, length 50m

PEG 20000 modified with 2-Nitroterephthalic acid

Max. Temperature: 220°C

Dimensionen

50 m x 0,32 mm ID

Artikel Nr.

[232084](#)

Probe

Alkohole (C1-C5)

Substanzen: Methanol / Ethanol / n-Propanol / Butanol-2 / n-Butanol / n-Pentanol


Ofen

80°C isotherm

Trärgas

0,55 bar N2 x Split 1:40

Detektor

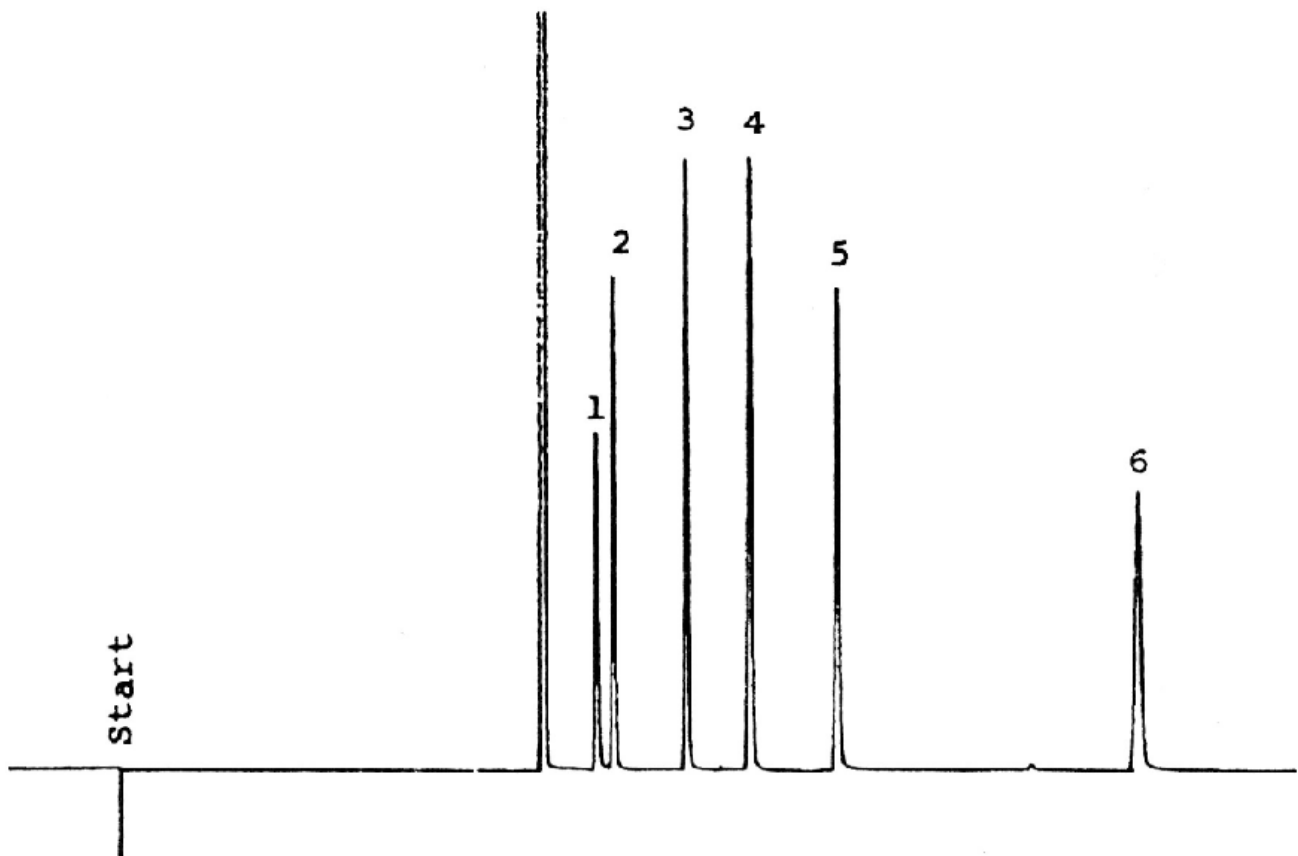

1.) Methanol

4.) Butanol-2

2.) Ethanol

5.) n-Butanol

3.) n-Propanol

6.) n-Pentanol