



Belegung

FS-OV-1

Dimensionen

25 m x 0,32 mm ID

Artikel Nr.

[232107](#)

Probe

Alkane/Aromaten

Substanzen: Pentan-isomere / n-Hexan / Benzol / Cyclohexan

Ofen


80°C isotherm

Trägergas

0,5 bar N<sub>2</sub> x Split 1:50

Detektor

FID 250°C, 10 x 16

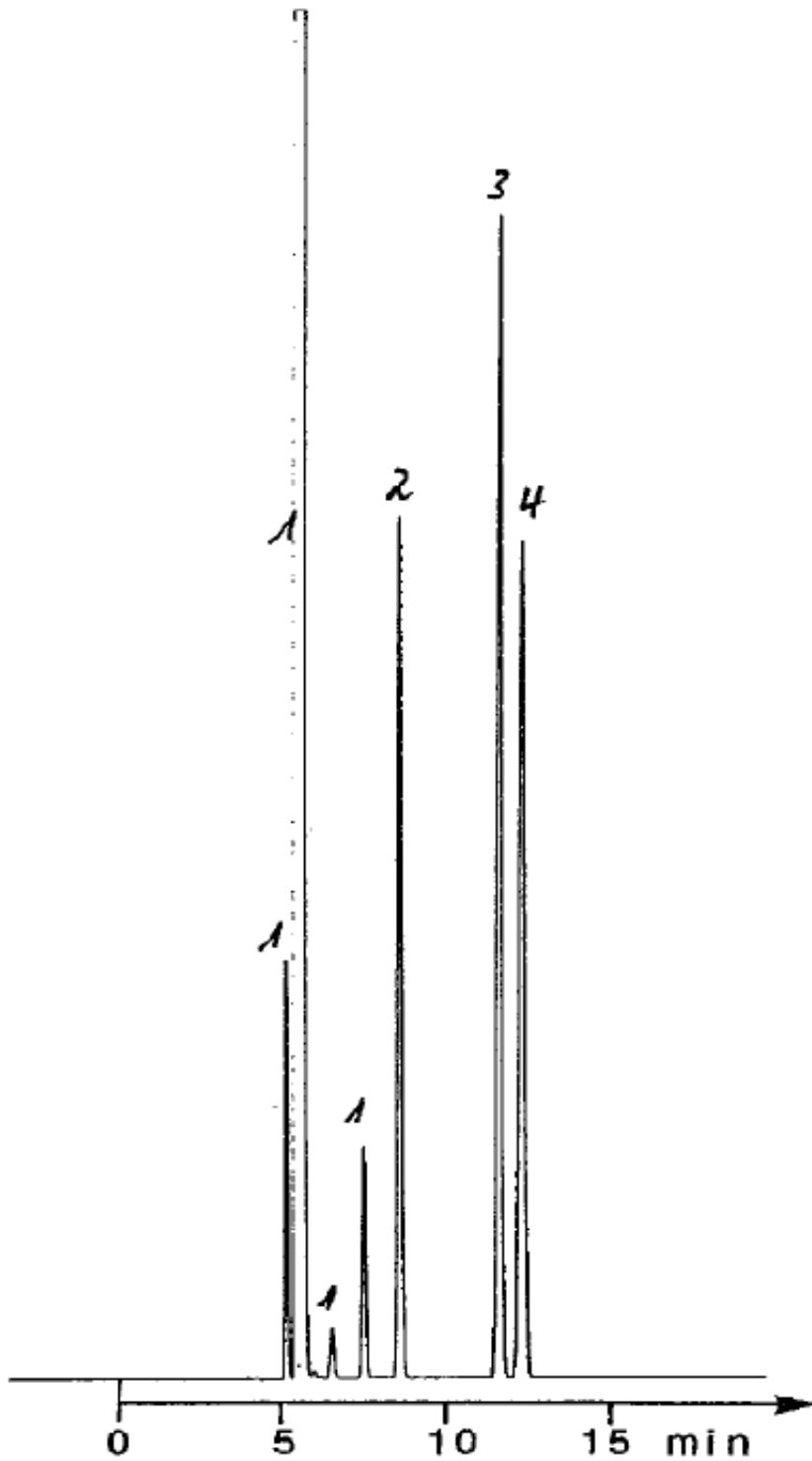

1.) Pentan-Isomere

2.) n-Hexan

3.) Benzol

4.) Cyclohexan