




Belegung	FS-SE-54 <i>ID 0,32mm, df 5µm, length 50m</i> <i>94% Methyl- 5% phenyl- 1% vinylsilicone</i> <i>Max. Temperature: 300/320°C</i>
Dimensionen	50 m x 0,32 mm ID
Artikel Nr.	<a href="#">232105</a>
Probe	Ethylenoxid und Acetaldehyd in Wasser
Ofen	80°C isotherm
Trägergas	1 bar N2 x Split 1:50
Detektor	FID 250°C, 10 x 8

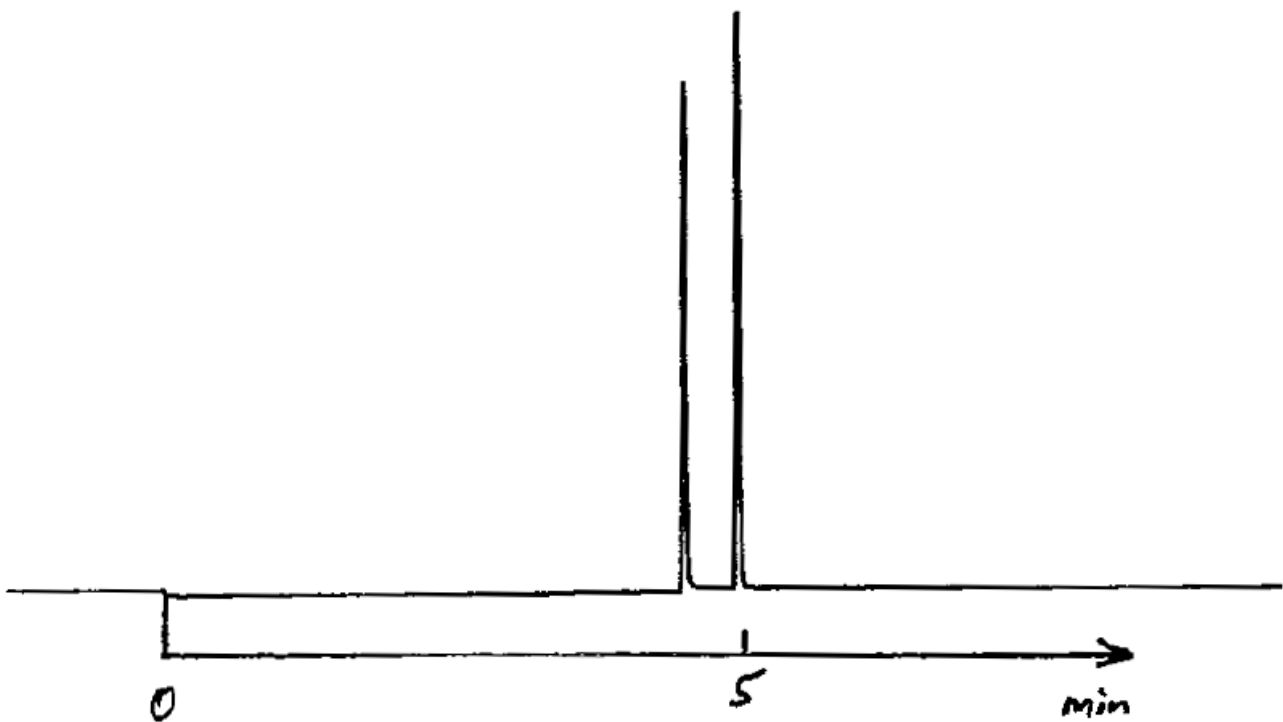

1.) Ethylenoxid

2.) Acetaldehyd