




Belegung	FS-SE-54 <i>ID 0,32mm, df 5µm, Länge 50m</i> <i>94% Methyl- 5% phenyl- 1% vinylsilikon</i> <i>Max. Temperatur: 300/320°C</i>
Dimensionen	50 m x 0,32 mm ID
Artikel Nr.	232105
Probe	Ethanol/iso-Propanol in Wasser
Ofen	80°C isotherm
Trägergas	1 bar N2 x Split 1:30
Detektor	FID 250°C, 10 x 32

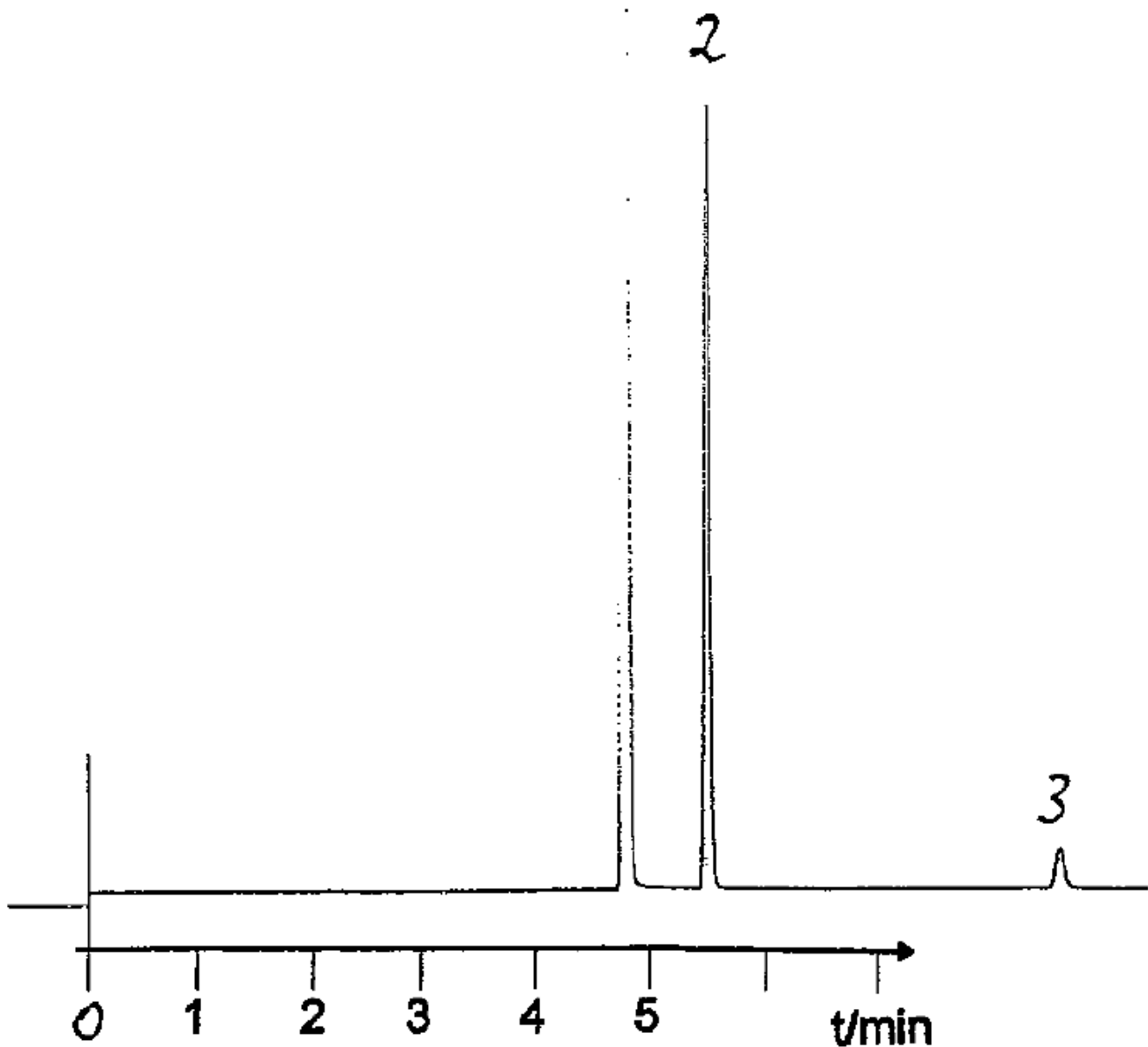

1.) Ethanol

2.) Isopropanol

