




Belegung	FS-OV-1 <i>ID 0,32mm, df 1µm, length 25m</i> <i>100% Dimethylsilicone, gum</i> <i>Max. Temperature: 300/320°C</i>
Dimensionen	25 m x 0,32 mm ID
Artikel Nr.	<a href="#">232023</a>
Probe	Essigsäure/Acetamid/Acetonitril in Wasser
Ofen	80°C isotherm
Trärgas	1 bar N2 x Split 1:30
Detektor	FID 250°C, 10 x 32

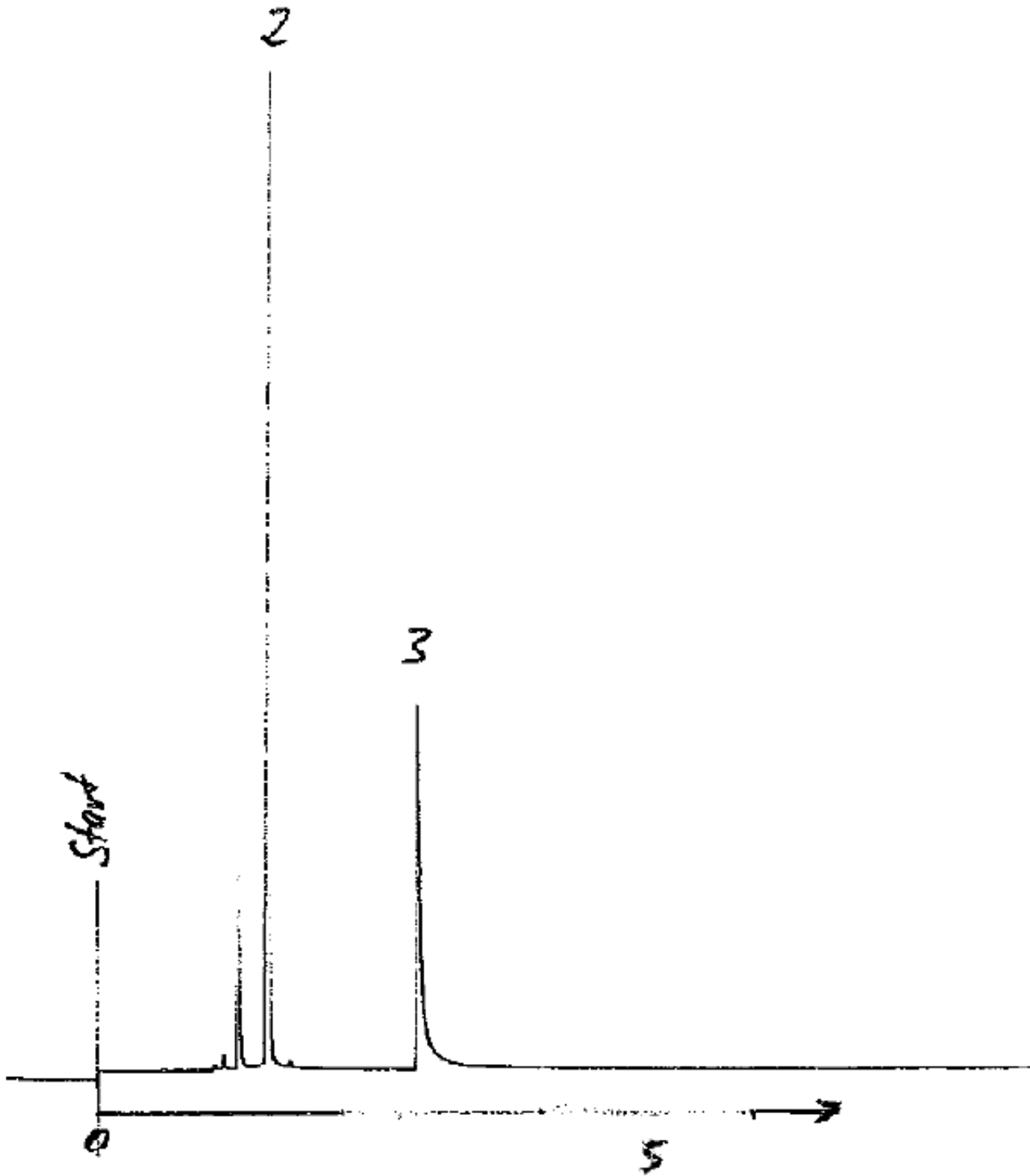

1.) Essigsäure

3.) Acetamid

2.) Acetonitril